

Single mode, up to 300mW, 1nm, Benchtop or Module



DATASHEET

Return to the Webpage



Specifications

Parameter Min **Typical** Unit Max Threshold Current 150 mΑ **Forward Current** Α Forward Voltage 2 ٧ **Optical Output Power** 300 mW/ Center Wavelength@25°C 1420 1510 nm Spectral Linewidth (FWHM) 1 nm Rise Time in Pulse Mode 30 ns Monitor Current @ VrPD=5V 1 mA PD reverse voltage 20 ٧ **TEC Current** 4 Α **TEC Voltage** 4.5 ٧ Operating Temperature (heatsink temperature) -5 +70 $^{\circ}\text{C}$ Storage Temperature -40 +75 °C

Agiltron provides cost-effective fiber-coupled laser sources with a wide range emitting spectrum from 370nm to 2000nm and line width from 10kHz to broadband to select. Each benchtop laser source features a pigtailed laser and high-precision, low-noise auto-feedback drive electronics to ensure constant output power or a constant driving current, and an integrated temperature control unit maintains optimal operating conditions. Each unit features a front fiber output connector and a universal power supply compatible with 100 to 240 VAC. We offer two packages: benchtop for ease of use and compact module for system integration. The user interface benchtop includes an intuitive LCD display for independent control of output power and temperature via two front rotating knobs. The module has two front output power and temperature settings. All units have a built-in isolator option to prevent reflection-induced laser emissions instability. We produce fiber-coupled isolators from 370nm to

2000nm. An isolator is essential to obtain stable laser output.

Features

- Turnkey Laser Source
- High Stability
- Advanced Feedback Control

Applications

- Medical Laser Treatment
- Biotechnology
- Others



Note: The specifications provided are for general applications with a cost-effective approach. If you need to narrow or expand the tolerance, coverage, limit, or qualifications, please [click this <u>link</u>]:

Rev 09/21/25

© Photonwares Corporation



Single mode, up to 300mW, 1nm, Benchtop or Module



DATASHEET

Benchtop Laser Source Operation Manual



- Plug in power cable
- Turn on Power Switch
- Setting the Output Power by rotating the knob
- Setting the laser diode Temperature by rotating the knob
- Connect a fiber path cable with matching connector (FC/APC is the default)
- Push the Emission switch to turn on the laser
- Measure the output power to verify

Module Laser Source Operation Manual



- Plug in power cable
- Turn on Power Switch
- Setting the Output Power by rotating the screw
- Setting the laser diode Temperature by rotating the screw
- Connect a fiber path cable with matching connector (FC/APC is the default)
- Push the Emission switch to turn on the laser
- Measure the output power to verify



Single mode, up to 300mW, 1nm, Benchtop or Module




DATASHEET

Module (mm): Plug-Play (Power Supply Included)






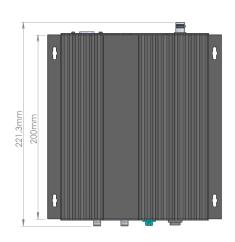



Benchtop (mm): Plug-Play 100-240VAC, USB

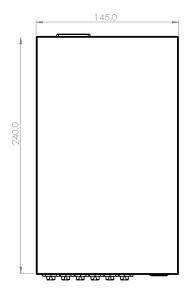





*Product dimensions may change without notice. This is sometimes required for non-standard specifications.

© Photonwares Corporation

P +1 781-935-1200

E sales@photonwares.com


w www.agiltron.com



Single mode, up to 300mW, 1nm, Benchtop or Module



Typical Spectrum

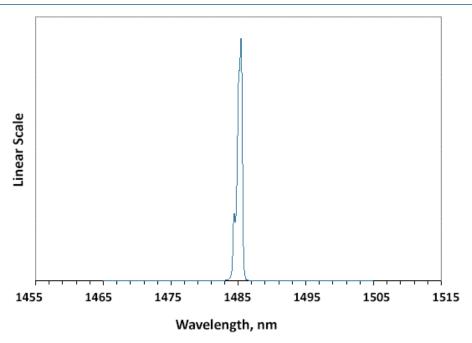

Ordering Information

Prefix	Wavelength	Power	Linewidth	Package	Isolator	Control Mode	TEC Cooling	Fiber Type	Connector
FCLS-	1485nm = 1485	180mW = 1 200mW = 2 220mW = A 240mW = B 260mW = C 280mW = D 300mW = 3	1nm = 1	Benchtop = 1 Module = 2	None = 1 Yes = 2	Constant Current = 2 Constant Power = 1	No = 1 Yes = 2	SM28 = 1 PM1550 = 5 50/125 = M Special = 0	FC/APC = 3 FC/PC = 2 Non = 1 SC/PC = 4 SC/APC = 5 LC/PC = 7 LC/UPC = U Special = 0